

PACKING INSTRUCTION




Important information about packaging and disposal


Don't disposed the
packaging after incoming.



Verify the packaging of
damages.

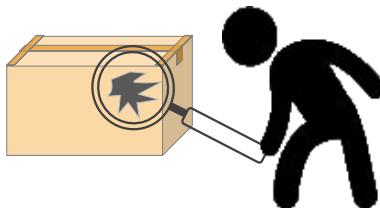

Retain the packaging in
case you need to ship the
device in the future.

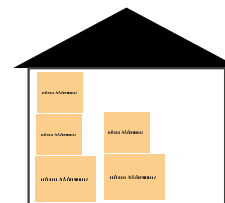



You will find the packing instruction on the page 2



Packing instruction

The packaging can be used for the devices below. Observe the instruction concerning the dimensions and weight of devices.

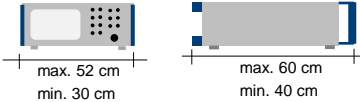




min. 13 cm
max. 18 cm

max. 52 cm
min. 30 cm

max. 60 cm
min. 40 cm

Weight of the devices from 15 kg up to 30 kg max.



CARDBOARD CUTTING 5001.7330.00
(at least height compensation 10 cm)

Check it out by pressing from above of the corrugated cardboard blanks whether it build up still a empty space above in the cardboard box.
If yes you have to fill up further corrugated cardboard blanks inside the cardboard box.

5001.7300.00	CARDBOARD MEM LARGE (1x)
5001.7317.00	CARDBOARD FOIL MEM LARGE. (2x)
5001.7323.00	CARDBOARD CUT MEM LARGE. (2x)
5001.7330.00	CARDBOARD CUT INSERT MEM LARGE
5027.8940.00	CARDBOARD HARDENING MEM LARGE



If you have some accessories like a power cable.
Pack it accordance the pictures below.



Apply the power cable into a bag between the membranes foils.
Fix them by adhesive tape on the surfaces of the membranes.



Apply the power cable into the membrane (empty space) above.

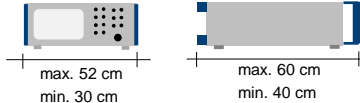




min. 17 cm
max. 22 cm

max. 52 cm
min. 30 cm

max. 60 cm
min. 40 cm

Weight of the devices up to 30 kg max.



CARDBOARD CUTTING 5001.7330.00
(at least height compensation 4 cm)

Check it out by pressing from above of the corrugated cardboard blanks whether it build up still a empty space above in the cardboard box.
If yes you have to fill up further corrugated cardboard blanks inside the cardboard box.



min. 22 cm
max. 25 cm

max. 52 cm
min. 30 cm

max. 60 cm
min. 40 cm

Weight of the devices up to 30 kg max.



Check it out by pressing from above of the corrugated cardboard membrane whether it build up still a empty space above in the cardboard box.
If yes you have to fill up further corrugated cardboard blanks inside the cardboard box.

Position of the devices in the membran packaging



The position of the devices into the membrane packaging has to packed centered with the handles to the short front side of the cardboard box.



The device into the membrane packaging is packed wrong.